
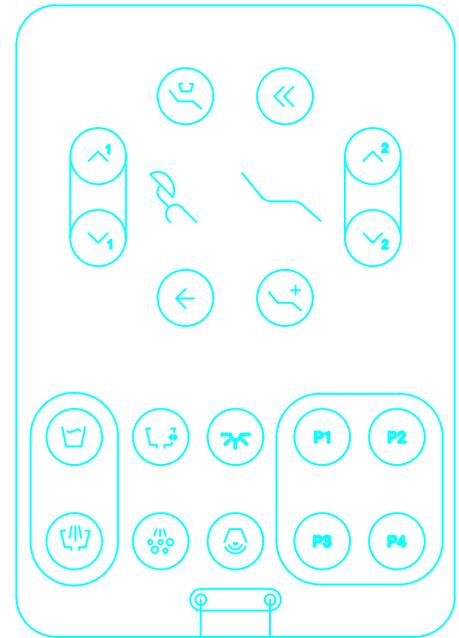


Wiring diagram of the Assistant Panel - Modul in assistant panel used for
MODEL PRO 500 (DC 310), MODEL PRO 700 (DC 350), DC 240, DC 290, DL 210, DL 270, DL 290, DL 320, DL 330, DL 340, DL 380,
MODEL PRO 600 (DA 270), DA 280, MODEL PRO 800 (DA 370), DA 380
 (11.5.2021)

- LEDs:**
 LD1 (GN) ~ Reserve
 LD2 (GN) ~ +3V3
 LD3 (GN) ~ System running indication
 LD4 (GN) ~ Reserve
 LD5 (GN) ~ Big Aspirator In
 LD6 (GN) ~ Reserve
 LD7 (GN) ~ Small Aspirator In
 LD8 (GN) ~ Syringe In
- JP1,3 ~ Temporarily unused (OFF)
- X1 ~ Not used
 X7 ~ Not used
- Turn ON JP2 -> Restart Dental Unit
 -> Wait for 5 min -> Turn OFF Dental Unit -> Turn OFF JP2 -> Turn ON Dental Unit)

Bowl Flushing Cup Filling Big Aspirator Small Aspirator

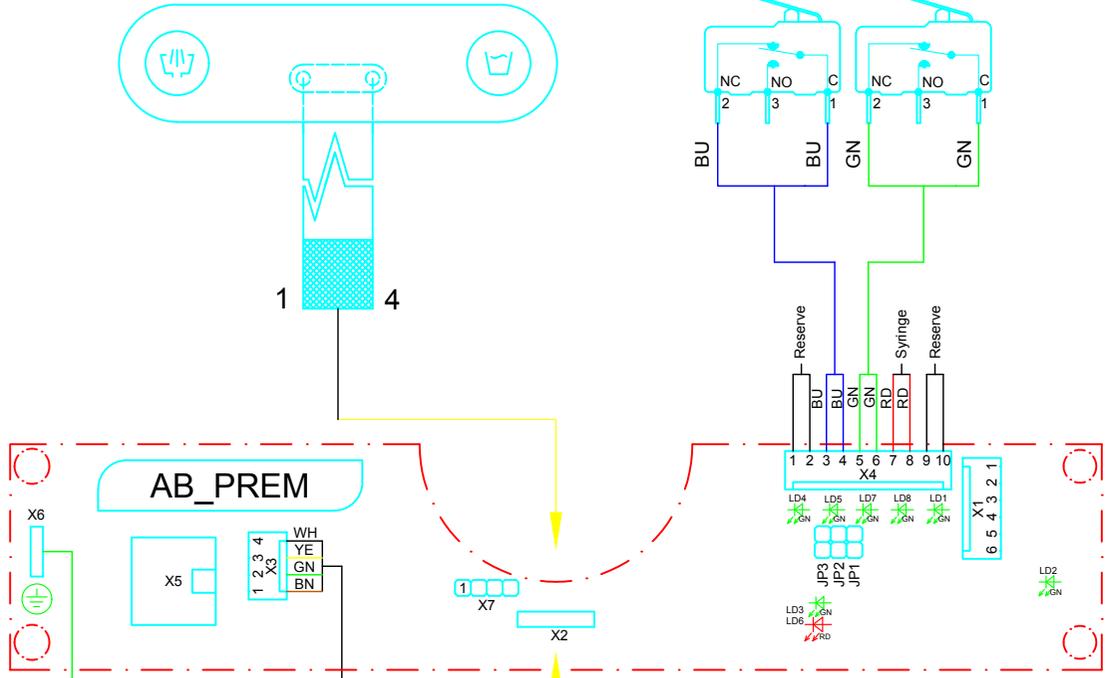

Assistant panel

Spittoon Block

SB_PREM
X9